

HOMECARE SERIES

BiliTx[™] **Biliblanket Care** *Training and Safety Guidelines*

For Your Baby's Safety

Following these safety guidelines will help reduce potential risks. Refer to the manufacturer website if you would like further information about your equipment.

- Use the biliblanket only for its intended use.
- Follow your physician's advice when using oxygen with phototherapy.
- Do not operate the biliblanket without the T-Vest (D).
- Do not spill liquid onto the equipment. Electric shock hazard and or equipment failure could result.
- Do not leave your baby unattended during phototherapy.
- Allow air to circulate around the unit.
- Turn the Illuminator Unit (A) off when the baby is not in the T-Vest (D).

When setting up the biliblanket, choose a place with the following features:


- Room temperature at 72°F or higher.
- A sturdy, flat surface away from drafts (bed, crib, or dining table).
- An accessible electrical outlet.
- Easy access for continual caregiver attendance.
- Absence of flammable substance or cleaning agents.



Setting up Your Biliblanket

Equipment Contents

- A. Illuminator Unit
- B. Ferrule
- C. Power cord
- D. T-Vest
- E. Panel

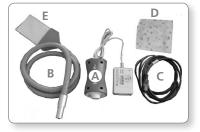

- 1 Place the **Illuminator Unit (A)** on a hard, flat surface. You may also hang the unit on the outside of a crib using the carrying case.
- 2 Insert the Ferrule (B) with the metal peg facing up, into the opening of the Illuminator Unit (A). Twist the Ferrule (B) one-quarter turn to lock into place.



- **3** Plug the **Power Cord (C)** into an electrical outlet.
- 4 Slide the **Panel** (**E**) with the sticker side down all the way into the **T-Vest** (**D**) pocket.
- 5 Lay the covered Panel (E) on a flat surface making sure the sticker side of the panel is facing down.





Placing your Baby in the T-Vest

1 Lay the baby on top of the T-Vest (D) covered Panel(E). The baby's back should be on the Panel (E) with the cable laying between the baby's legs.



2 Wrap the T-Vest (D) around the baby and secure with tab.



T-Vest Tips

- Make sure the T-Vest (D) is not wrapped too tight. Your finger should easily fit between the T-Vest (D) and the baby's body.
- Wrap your baby in a blanket or put the baby in a sleeper for a tighter fit.
- For a larger or more active baby, tape the T-Vest (D) to the baby's diaper.
- If the T-Vest (D) gets soiled, discard it and use a new one.

Starting Therapy

- **1** Turn the Illuminator Unit (A) on by pressing the ON button.
- **2** Make sure the green power light is on.
- **3** Verify the Panel (E) is lighting under the baby.

Call Intermountain Homecare & Hospice if...

• You experience any problems with the equipment.

Troubleshooting

Problem	Action
Therapy ON/OFF button is not green	 Make sure power cord is properly attached and plugged into an active electrical outlet Verify the unit is turned on
Fiber optic panel light is not working	Make sure the panel is securely locked into illuminator device
Yellow LED is flashing	Make sure the panel is securely locked into the illuminator
Device will not turn off with the ON/OFF button	Press and hold the Therapy ON/OFF button for 3 seconds



Content created and approved by Intermountain Homecare & Hospice