


Giving an Enteral Feeding with an Infinity Pump


Priming and loading the disposable tubing set:

- 1 Prior to setting up the pump, fill the feeding bag with formula, gently squeeze the air out of the bag, and close the cap.

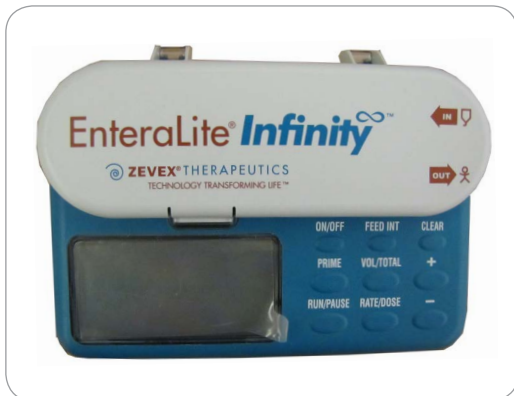

- 2 Remove the protective cover from the end of the tubing.

Use the white cone adapter or remove it and use the **ENFit** connection, as needed, to connect to your tube.

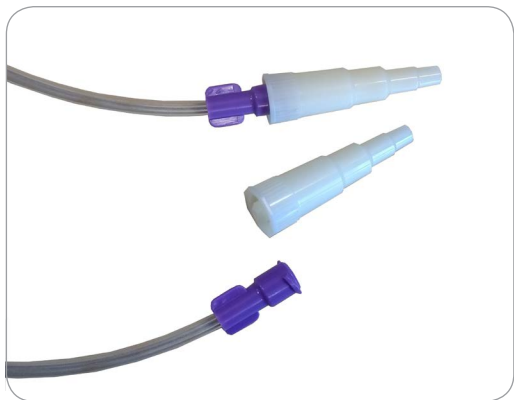

- 3 Gently pinch and hold the teal-colored tubing below the teardrop symbol. This will open the valve and let you squeeze the tube with your hands to allow the formula in.



- 4 Open the pump door and loop the teal-colored silicone tubing around the pump wheel, stretching it lightly. Seat the cassette into the pump.


Close the pump door. The door must be closed before the pump will work. If the door breaks, call the Enteral Team.



Attach the pump to the IV pole and hang the feeding bag.

OR

Fit the pump, bag and tubing into a backpack, securing the pump, bag, and tubing with the straps.

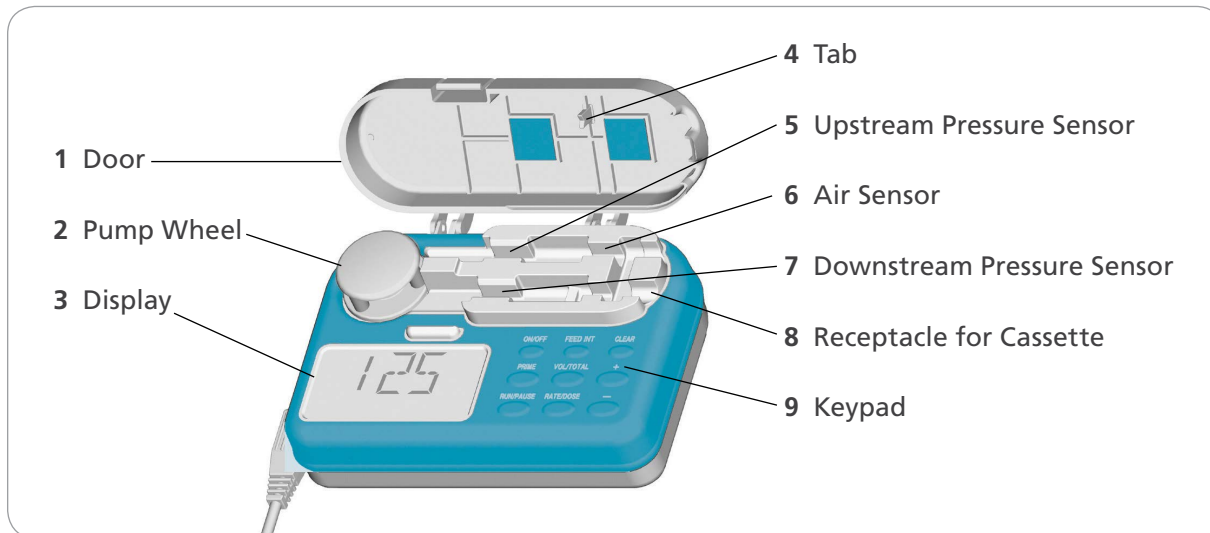


Programming and running the pump:

1 Press and hold the **ON/OFF** key for 1.5 seconds to turn the pump on. Numbers will flash across the screen — this is an automatic self-test. After running the self-test, the display will show the last programmed rate.




2 Press the **RATE/DOSE** key to display the rate (see the bottom left corner of the screen). Program the prescribed rate using the **(+)** and **(-)** keys.

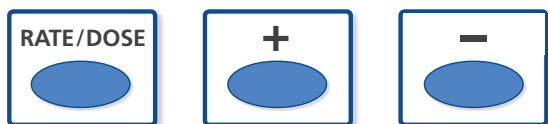

3 Press the **RATE/DOSE** key to display the dose.

Program the dose using the **(+)** and **(-)** keys. Set the dose to **INF**. This means the pump will not alarm **DOSE DONE**. The pump will run continuously as long as the feeding solution remains in the bag. When empty, the pump will alarm **NO FOOD** or **NO FLOW IN**.

Flush the feeding tube.

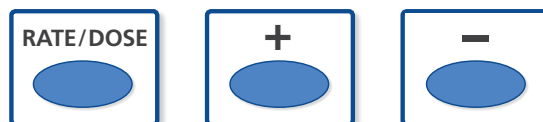

4 Attach the primed feeding bag and tubing to the feeding tube.

5 Press the **RUN/PAUSE** key to start the pump.



6 When feeding is complete, press the **ON/OFF** key for 1.5 seconds to turn the pump off.



7 Flush your feeding tube after disconnecting.

See troubleshooting tips on page 3.

Troubleshooting guide for alarms

<p>No Flow In OR No Flow Out</p>	<p>No Flow In (blocked between delivery set and pump) OR No Flow Out (blocked between pump and patient)</p> <ul style="list-style-type: none"> • Is delivery set tubing pinched, kinked or clogged? Remove kinks, pinches, or clogs. Resume feeding. • If the pressure sensor area is soiled, clean it with a cotton swab, soft cloth, dampened sponge, or wash the pump under running water. • A broken tab on the door may cause this alarm. Call the Enteral Team.
<p>ERRA, ERRD, ERRF</p>	<ul style="list-style-type: none"> • May happen if the pump door is open or was opened while the pump was on. Turn the pump off, open the door, unload and reload cassette, close the door, and then turn the pump back on.
<p>Load Set</p>	<ul style="list-style-type: none"> • Turn the pump off. • Check that the door is closed. • Unload and reload the cassette, then close door and turn the pump back on.
<p>No food alarm</p>	<ul style="list-style-type: none"> • Rinse and refill the bag with formula if it is empty. • If air is in the bag or the tubing, squeeze it out of the tube with your hand and resume the feeding. • If the formula is foamy or has bubbles, let it sit for 10 to 15 minutes before proceeding with the feeding.
<p>Low battery</p>	<ul style="list-style-type: none"> • Is the A/C adapter charger plugged into a wall outlet and into the pump properly? Make sure the outlet is functioning properly. The green light on the power cord will light up and the battery symbol on the pump will flash up and down.
<p>Charging problems</p>	<ul style="list-style-type: none"> • Check the connection for both the wall plug and pump adapter. A green indicator light on the charger wall plug should be lit and the pump display should have a plug symbol with scrolling bars between E and F of the fuel gauge. When the pump is fully charged, the bars will stop scrolling. • Plug another device into the outlet to verify that the outlet is working properly. • If a battery symbol and E and F of the fuel gauge is flashing, contact the Enteral Team.

Need help or have questions?

Enteral Team
385-887-7356

Intermountain Homecare Main Numbers
385-887-6000 | 800-527-1118 (Toll Free)

Content created and approved by Intermountain Homecare & Hospice

Intermountain Healthcare complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. Se proveen servicios de interpretación gratis. Hable con un empleado para solicitarlo. 我們將根據您的需求提供免費的口譯服務。請找尋工作人員協助。