

BiliTx™ Biliblanket Care: Training and safety guidelines

For your baby's safety

Following these safety guidelines will help reduce potential risks.

- Use the biliblanket only for its intended use.
- Follow your doctor's advice when using oxygen with the biliblanket.
- Do not operate the biliblanket without the T-vest (D).
- Do not spill liquid onto the equipment. Electric shock and/or equipment failure could result. If it does get wet, unplug the biliblanket immediately and dry off. Turn the biliblanket back on when it is fully dry.
- Do not leave your baby unattended while using the biliblanket.
- Allow air to circulate around the illuminator unit.
- Turn the illuminator unit (A) off when the baby is not in the T-vest (D).

Do not use any lotions, oils, or creams on the baby's exposed skin during therapy.

Refer to the manufacturer's website if you would like further information about your equipment.

When setting up the biliblanket, choose a place with the following features:

- Room temperature at 72°F or higher.
- A sturdy, flat surface away from drafts, such as a bed, crib, or dining table.
- Easy access to an electrical outlet.
- Easy access for caregivers.
- Away from flammable substances or cleaning supplies.



Setting up your biliblanket

The biliblanket kit includes a(n):


A. Illuminator unit


B. Ferrule

C. Power cord

D. T-vest


E. Panel

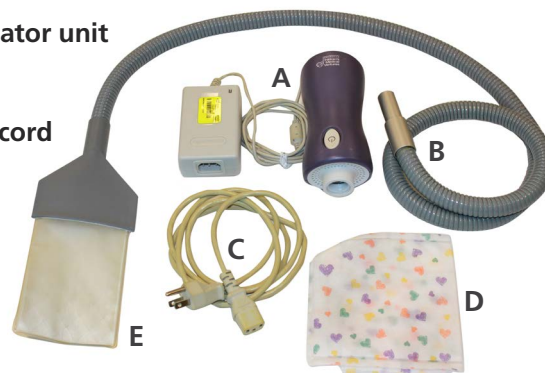

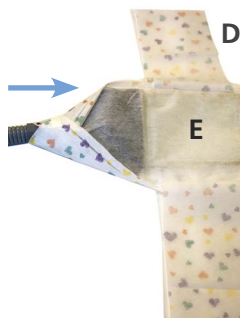- 1 Place the **illuminator unit (A)** on a hard, flat surface. You may also hang the unit on the outside of a crib using the carrying case.
- 2 Insert the **ferrule (B)** with the metal peg facing up, into the opening of the **illuminator unit (A)**. Twist the **ferrule (B)** one-quarter turn to lock into place.
- 3 Plug the **power cord (C)** into the **illuminator unit (A)** and an electrical outlet.
- 4 Slide the **panel (E)** with the sticker side down all the way into the **T-vest (D)** pocket. 
- 5 Lay the covered **panel (E)** on a flat surface making sure the sticker side of the panel is facing down.

Placing your baby in the T-vest

1 Lay your baby on top of the **T-vest (D)** covering the **panel (E)**. Keep your baby's diaper on. Your baby's back should be on the **panel (E)** with the cable laying between your baby's legs.

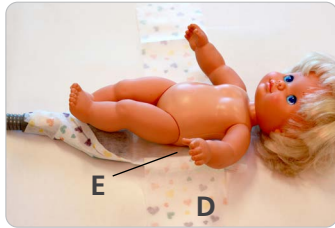

2 Wrap the **T-vest (D)** around your baby (one side at a time) and secure with the tab.



T-vest tips

- Make sure the **T-vest (D)** is not wrapped too tight. Your finger should easily fit between the **T-vest (D)** and your baby's body.
- Wrap your baby in a blanket or put your baby in a sleeper for a tighter fit. During therapy, you may lay your baby in their crib or hold them in your arms.
- For a larger or more active baby, tape the **T-vest (D)** to your baby's diaper.
- If the **T-vest (D)** gets soiled, discard it and use a new one. You will receive a few extra **T-vests (D)**.

Starting therapy

- 1** Turn the **illuminator unit (A)** on by pressing the ON button.
- 2** Make sure the green power light is on.
- 3** Verify the **panel (E)** is lighting under the baby.

**Call Intermountain
Homecare & Hospice
if you experience any
problems with the biliblanket:
1-800-527-1118**

Troubleshooting

Problem	Action
Therapy ON/OFF button is not green	Make sure the power cord is properly attached to the illuminator unit and plugged into an active electrical outlet. Double check that the unit is turned on.
Panel is not lighting	Make sure the ferrule is securely locked into the illuminator unit.
Yellow LED light on the illuminator unit is flashing	Make sure the ferrule is securely locked into the illuminator unit.
Device will not turn off with the ON/OFF button	Press and hold the Therapy ON/OFF button (located on the illuminator unit) for 3 seconds.

Content created and approved by Intermountain Homecare & Hospice

Intermountain Healthcare complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. Se proveen servicios de interpretación gratis. Hable con un empleado para solicitarlo. 我們將根據您的需求提供免費的口譯服務。請找尋工作人員協助。

© 2010-2021 Intermountain Healthcare. All rights reserved. The content presented here is for your information only. It is not a substitute for professional medical advice, and it should not be used to diagnose or treat a health problem or disease. Please consult your healthcare provider if you have any questions or concerns. Intermountain Patient Education FSHC005 - 11/21 (Last reviewed - 11/21) Also available in Spanish.