

Giving Your Medicine with a Mini-Bag Plus

Your doctor ordered a medicine that you will give in your intravenous (IV) line. This medicine comes in a Mini-Bag Plus, which is a vial of medicine connected to a small bag of sterile IV solution. You will mix the vial into the bag and then you give the medicine in your IV.

You will also need to “flush” your IV line to keep it clear. Flush means you fill the IV tubing with solution to keep your IV line from clotting. Your nurse will teach you how to prepare and give your medicine and flush your line safely.

Follow the steps in this fact sheet and your nurse’s instructions carefully. Unless your nurse tells you otherwise, discard all used supplies in the regular trash.

Medication information

Medication: _____

Schedule: _____

S.A.S. FLUSH

Saline: _____ mL

Administer medication

Saline: _____ mL

S.A.S.H. FLUSH

Saline: _____ mL

Administer medication

Saline: _____ mL


Herapin: _____

Prepare

- 1 | **Clean** your hands with cleansing solution.
 - **Rub** your hands vigorously until they are dry.
 - **Clean** your hands each time you work with your IV.



- 2 | **Check** the medication label.
 - **Verify** each label has your name and dose, and is not expired.

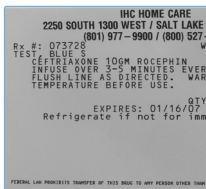

- 3 | **Lay** out your supplies.
 - Mini-Bag with medicine
 - Warm to room temperature if refrigerated.
 - If medicine in the vial is not dry, call the pharmacist.
 - IV Tubing
 - Alcohol disinfectant cap
 - Sterile wipes
 - Flush syringes
 - 2 Saline syringes
 - Heparin syringe, if needed



Prepare your dose

- 1 | **Bend** the tube connecting vial to bag, up and then down to break seal.
- 2 | **Fill** the vial by holding the bag with the vial down and squeezing the bag until the vial is half full.

Swirl gently to mix.

- 3 | **Turn** the bag so the vial is on top and squeeze bag to force air into vial — this empties the vial into the bag.

- 4 | All but a little of the medicine in the vial should empty into the bag.


Repeat if needed.

Leave the vial connected to the bag.



Connect IV tubing to the medicine bag


- 1 | **Close** the tubing clamp.

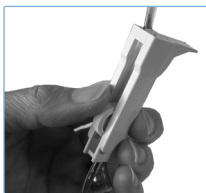


- 2 | **Pull** the cover off the medicine bag.



- 3 | **Remove** the spike cap.

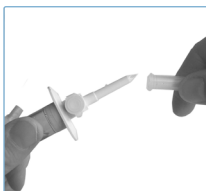


- 4 | **Push** the spike into the medicine bag.

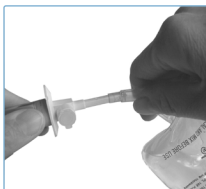

Prime your IV tubing

- 1 | **Hang** the bag on the IV pole.

- If there is no fluid in the drip chamber, **gently squeeze** and release it to add a little fluid.

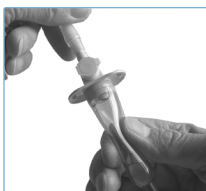

- 2 | **Open** the IV tubing clamp.



- 3 | **Fill** tubing with the solution.

- **Loosen** the tubing cap if solution is not flowing, then retighten.



- 4 | **Close** the clamp.

- **Reclamp** the tubing, if necessary.




Flush your IV line


Flush following the “How to flush your IV line” instructions.

How to flush your IV line

- 1 | **Remove** the alcohol disinfectant cap and discard.

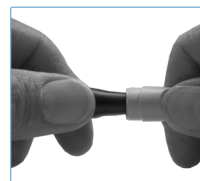

- 2 | **Attach** the syringe firmly to the IV needleless connector using a push and twist motion.

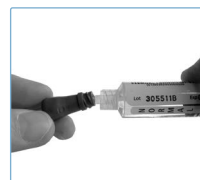


- 3 | **Open** the clamp.




- 4 | **Inject** the solution.




- 5 | **Close** the clamp.

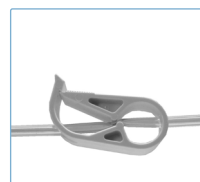

- 6 | **Remove** and discard the syringe.

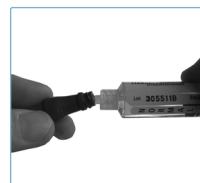

- 7 | **Attach** a new alcohol disinfectant cap after the final flush.

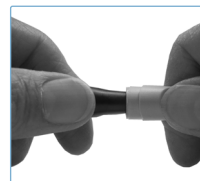


* If the needleless connector touches anything, **CLEAN IT** with a sterile wipe for 15 seconds.

Administer your medication

- 1 | **Scrub** the needleless connector with a sterile wipe for 15 seconds if it has touched anything.



- 2 | **Attach** the medicine bag firmly to the needleless connector.

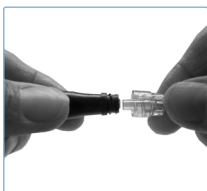

- 3 | **Open** the IV line clamp.



- 4 | **Adjust** the rate to _____.

- 5 | **Infusion** will take _____ minutes.

- 6 | **Stop** your infusion before the bag is completely empty.


After your infusion

When your infusion is done, a small amount of fluid will be left in the bag. This is considered overfill and doesn't affect your dose. Follow the steps below once your infusion is complete.

- 1 | **Clean** your hands with cleansing solution.



- 2 | **Disconnect** IV tubing from your IV line.

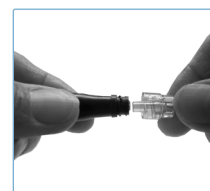

Reuse the tubing


Tubing may be reused for up to 24 hours. Follow the steps below to reuse your tubing. Discard your tubing after 24 hours.

- 3 | **Protect** the end of the tubing with a new alcohol cap.



- 4 | **Leave** the tubing attached to the empty bag until the next dose.

- With your next infusion, **remove** old bag and **connect** new bag to your tubing.
- **Reprime** the tubing if necessary.

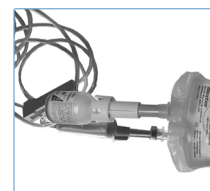

Need help or have questions?

Nurse: _____

Phone number: _____

Intermountain Healthcare & Hospice
385-887-6000
800-527-118

Flush Your IV Line

Follow the "How to flush your IV line" instructions on page 2:

- **Flush** with saline.
- **Flush** with Heparin, if instructed.

Notes

[illegible]

Intermountain Health complies with applicable federal civil rights laws and does not discriminate on the basis of race, color, national origin, age, disability, or sex. Se proveen servicios de interpretación gratis. Hable con un empleado para solicitarlo. 我們將根據您的需求提供免費的口譯服務。請找尋工作人員協助。

© 2014-2024 Intermountain Health. All rights reserved. The content presented here is for your information only. It is not a substitute for professional medical advice, and it should not be used to diagnose or treat a health problem or disease. Please consult your healthcare provider if you have any questions or concerns. FSHC023 - 12/24 (Last reviewed - 02/23) Also available in Spanish.